





### HUMIDITY SEPARATORS **S11** (Centrifugal design)

### DESCRIPTION

When wet steam is used in sterilization, moisture in suspension reduces the heat transfer efficiency, and the validity of the sterilization process can be compromised.

The ADCAPure S11 series centrifugal separators remove moisture from steam pipelines. Steam passes through the separator and, as a result of centrifugal forces, impact and swirling effects, the particles with a heavier specific gravity are separated, such as water droplets and moisture in suspension.

The condensate collected at the bottom of the separator must be automatically drained by a suitable steam trap.

#### MAIN FEATURES

316L stainless steel construction. No moving parts. Self drainable design. Clamped construction to ease access for cleaning operations.

#### STANDARD SURFACE FINISH

Internal wetted parts: ≤ 0,51 µm Ra – SF1. External: Satin bead blast finish - 1,6 µm Ra. Other surface conditions see TIS.GIA - General information ADCAPure.

OPTIONS:	Vent connection. Special connections and dimensions.					
USE:	Steam, compressed air and other gases.					
AVAILABLE MODELS:	S11 – horizontal connections. S11A – S11 with air vent connection.					
SIZES:	1/2" to 2".					
CONNECTIONS:	ASME BPE clamp ferrules. Others on request.					
PACKAGING:	Assembling and packaging in a clean room certified according to ISO 14644-1. The product is end capped and sealed with recyclable thermo-shrinkable plastic film, to avoid contamination.					
INSTALLATION:	Always with the condensate discharge pointing downwards. See IMI – Installation and maintenance					

instructions.





CE MARKING – GROUP 2 (PED – European Directive)					
PN 10	Category				
1/2" to 1"	SEP				
11/2" to 2"	1				

<b>BODY LIMITING CONDITIONS *</b>					
ALLOWABLE PRESSURE	RELATED TEMPERATURE				
10 bar	50 °C				
8 bar **	175 °C				
7,4 bar	200 °C				


\* Other limits on request. Maximum operating conditions may be limited by the equipment end connections due to normative restrictions.

\*\* PMO – Max. operating pressure for saturated steam. Minimum operating temperature: -10 °C. Design code: AD-Merkblatt.

# VALSTEAM JDCJ

We reserve the right to change the design and material of this product without notice.



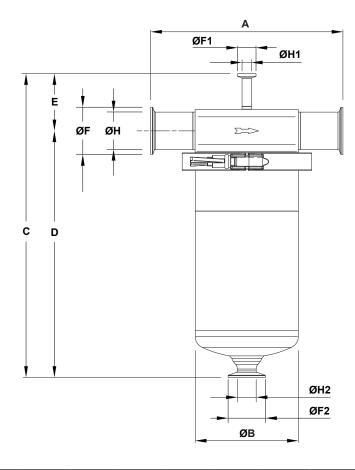

DIMENSIONS (mm)													
SIZE	Α	ØB	с	D	E	ØF	ØF1	ØF2	ØН	ØH1	ØH2	VOL. (dm³)	WEIGHT (kg)
1/2"	215	114	345	283,5	62,5	25	25	50,5	9,4	9,4	22,1	2,84	3,8
3/4"	215	114	345	283,5	62,5	25	25	50,5	15,8	9,4	22,1	2,87	3,9
1"	215	114	345	283,5	62,5	50,4	25	50,5	22,1	9,4	22,1	2,9	4,2
11/2"	235	141	416	338,5	77,5	50,4	25	50,5	34,8	9,4	22,1	5,82	7,25
2"	260	141	416	338,5	77,5	63,9	25	50,5	47,5	9,4	22,1	5,93	7,28


MATERIALS					
POS. Nº	DESIGNATION	MATERIAL			
1	Body	AISI 316L / 1.4404			
2	Cover	AISI 316L / 1.4404			
3	Clamp	AISI 316 / 1.4401			
4	* Gasket	** PTFE/FPM Envelope			
5	Internals	AISI 316L / 1.4404			

\* Available spare parts. \*\* Others on request.

EN 10204 3.1 certificate available if requested along with the order.

Remarks: FDA / USP Class VI seals certificate on request.

All separators have a serial number. In case of non-standard separators, this number must be supplied if spare parts are ordered.

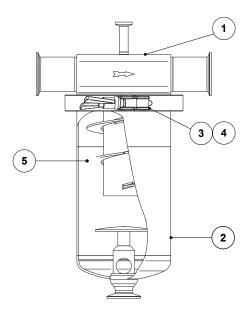

# VALSTEAM ADCA